

INSTRUCTION FOR USE

This device is indicated for specialized procedures, which must be carried out by professionals qualified in Implantology. For best results, use the product with appropriate techniques. Always apply the product under appropriate conditions, in a surgical environment.

INDICATIONS FOR USE

The Compact Surgical Kit consists of instruments that are used to apply specialized surgical procedures by qualified dental surgeons, including diagnosis, pre-operative planning and surgical protocol. The Compact Surgical Kit is indicated to assist in the installation of all different models of implants from DSP Biomedical.

COMPOSITION



Descrição	Código
Cylindrical Drill ϕ 1.6	28.8016C
Cylindrical Drill ϕ 2.0	28.8020C
Cylindrical Drill ϕ 2.5	28.8025C
Cylindrical Drill ϕ 2.8	28.8028C
Cylindrical Drill ϕ 3.0	28.8030C
Cylindrical Drill ϕ 3.15	28.8031C
Cylindrical Drill ϕ 3.3	28.8033C
Cylindrical Drill ϕ 3.5	28.8035C
Cylindrical Drill ϕ 3.8	28.8038C
Cylindrical Drill ϕ 4.3	28.8043C
Countersink Drill ϕ 4.0	28.8340
Countersink Drill ϕ 5.0	28.8350
Conical Drill ϕ 3.5	27.3515
Conical Drill ϕ 3.8	27.3815
Conical Drill ϕ 4.3	27.4315
Pilot Drill ϕ 2.0/ ϕ 3.0	28.8223C
Pilot Drill ϕ 3.0/ ϕ 3.8	28.8234C
Digital Key 0.9	28.2509

Digital Key 1.2	28.2512
Digital Key 1.3	28.2613
Ratchet Connection	28.4006
Ratchet/Contra Angle Connection	28.4008
Ti Hex Ratchet Connection (short)	28.2020
Ti Hex Ratchet Connection (long)	28.3020
Ti Hex Contra-Angle Connection	28.5020
Implants Meter	28.4015
Parallelizer	28.4041
Titanium Forceps	27.4040
Drill extender	28.4030
Surgical Torquimeter	28.4103
Autoclavable case	28.4105

APPLICATION

The Compact Surgical Kit consists of drills, installation wrenches and other accessories for the following purposes:

- Cylindrical Drill: used to deepen and direct the perforation of bone tissue, preparing the bone bed for implant installation. It has a presentation with Stop adaptation, the Stop does not come with the drills, it must be purchased separately. The Stop is an accessory that helps limit the depth of bone drilling.
- Countersink Drill: They are used as countersink drills at the end of the bone drilling sequence. Its purpose is to create a bevel for correct seating of the implant platform.
- Conical Drill: perforation of bone tissue during the preparation of the surgical socket prior to installation in bone types I and II (bone quality classification according to Lekholm and Zarb, 1985), they are used to progressively enlarge the diameter of the perforation until the dimension final is suitable for the installation of the conical implant. It has a presentation with Stop adaptation, the Stop does not come with the drills, it must be purchased separately. The Stop is an accessory that helps limit the depth of bone drilling.
- Pilot Drill: Used in the surgical sequence after the screw cutter and before the conical drill to mitigate bone heating.
- Digital Key: Recommended for the initial manual placement of prosthetic components.
- Ratchet Connection: Must be coupled to a Surgical or Prosthetic Torque Meter, it is recommended to finalize the placement/removal of a prosthetic component, in order to achieve a more precise torque.
- Ratchet/AC Connection: it is a connection that changes from ratchet to contra angle.
- Implant meter: it is a millimeter ruler where it is possible to measure the height of the implant part and the diameter of the drills.
- Parallelizer: Before use, it is recommended that a thread be passed through the central hole of the

parallelizer for handling inside the oral cavity. It is indicated to evaluate the angle at which bone drilling is being performed.

- Titanium forceps: used to pick up and transport implants.
- Drill extender: Also called extender, it is indicated to extend drills during surgical instrumentation, when the adjacent interdental height is greater than 20.0mm.
- Torquimeter: Used to torque the implant and/or prosthetic component during its insertion.
- Autoclavable case: Case for storing instruments for carrying out surgical procedures. It is made from autoclavable polymer. It has silicone supports for storing and securely fixing each instrument. It has markings that guide the sequence of use of the instruments and can be equipped according to the procedure. Product indicated for supporting and storing surgical and/or prosthetic instruments during their use and sterilization.

NOTICE

Failure to recognize actual drill lengths in relation to radiographic measurements can result in permanent damage to nerves and other vital structures. Drilling beyond the intended depth for lower jaw surgery may result in permanent numbness in the lower lip and chin or lead to bleeding in the lower part of the mouth. How to follow the mandatory procedures for any surgery, such as: asepsis during bone drilling, avoiding damage to blood vessels and nerves, using anatomical knowledge and pre-operative radiographs.

CONTRAINDICATIONS

This product has no contraindications, as long as it is used correctly for the indicated purposes.

RISKS AND BENEFITS

Like any surgery, it cannot be guaranteed one hundred percent, as achieving good performance involves several factors, including usability, the patient's clinical conditions and the product itself. Failure to observe the indicated usage limitations and work steps may result in failure.

Failure to recognize actual drill lengths in relation to radiographic measurements can result in permanent damage to nerves and other vital structures. Drilling beyond the intended depth for lower jaw surgery may result in permanent numbness in the lower lip and chin or lead to bleeding in the lower part of the mouth.

WARNINGS AND PRECAUTIONS

- Inadequate planning can compromise the performance of the implant/prosthesis set, resulting in system failures, such as loss or fracture of the implant, loosening or fracture of components and/or prosthetic screws.
- Due to its function, the drilling length must be a maximum of 0.5 mm greater than the implant insertion

depth. This additional length must be considered during the planning phase. • Drills must not be resharpened.

- Failure to replace the drills as recommended by the manufacturer may cause undue bone heating, compromising the success of the procedure.
- Care must be taken in cases of patients who show signs of allergy or hypersensitivity to the chemical components of the material: surgical stainless steel.
- Make sure to use the drill compatible with the sequence of drills indicated according to the dimensions and prosthetic interface of the planned implant.
- DSP Biomedical Drills are only compatible for preparation prior to the installation of DSP Biomedical implants.
- This product must be used sterile.
- To carry out the procedure, make sure that the patient has sufficient interocclusal space to handle the instruments in the desired region.
- Do not stop the motor rotation with the drill inside the surgical cavity, as this may make removal difficult or cause the drill to fracture.
- Do not use the product if its packaging is damaged.
- Do not use the product if its validity has expired.
- Before each procedure, check the perfect fit between the pieces.
- Ensure that the parts are not swallowed or aspirated by the patient.
- Make sure you have all the necessary instruments to perform the procedure according to the surgical plan.
- Before each procedure, check the condition of the instruments, always respecting their useful life. Replace instruments in case of damage, faded markings, compromised sharpening, deformation or wear.
- Always use the DSP Biomedical product sequence. The use of instruments and/or prosthetic components from other manufacturers does not guarantee the perfect function of the DSP Biomedical Implant System and exempts any product warranty.
- It is the dentist's responsibility to use DSP Biomedical products in accordance with the instructions for use.

IMPORTANCE OF THE NEED TO ADHERENCE TO A CARE REGIME

Products must be protected against aspiration when handled intraorally. Vacuuming products can lead to infection or unplanned physical injury. If you want to protect it, use rubber dam. If an implant or instrument has been swallowed or aspirated, call a doctor immediately. In addition to the mandatory precautions for any surgery, such as asepsis, during drilling into the jaw bone, damage to nerves and vessels must be avoided by referring to anatomical knowledge and pre-operative medical images (e.g. x-rays). Failure to recognize the

actual duration of drills in relation to radiographic measurements can result in permanent damage to nerves and other vital structures. Drilling beyond the depth intended for lower jaw surgery can potentially result in permanent numbness in the lower lip and chin or lead to hemorrhage in the floor of the mouth. Inappropriate use of products leads to poorly executed work and increased risk. In particular, users of hand tools should be careful to use them gently and with consideration. The user must always avoid touching instruments and parts without protection (sterile protective gloves and aprons must be used). Thermal bone damage caused by rotating and oscillating tools must always be avoided (user training, work at low speed and with sufficient cooling (see section "Cooling"). During intraoral application, attention must be paid to the fact that the products are protected from being sucked in or dropped to the ground. Rotary instruments need to be secured as far as possible with their set speed before applying them to the object are used with the rotary instruments. bone or fracture of system components. Drills are supplied in non-sterile conditions and must be reprocessed and sterilized before first use on the patient and immediately after each use, all products must be disinfected and sterilized. Inadequate cleaning and sterilization of instruments may result in patient infection with harmful bacteria. To avoid damaging instruments, they must be removed from the blister packaging individually. The shipping blister is not intended to be used as a container for steam sterilization of drills. They must be unpacked before the first reprocessing. It is essential to only use turbines, as well as manual and angular parts that are technically and hygienically impeccable, maintained and clean. Do not use the device if the primary package has been previously damaged or opened. Do not use damaged or blunt instruments for drilling. Broken instrument cutting edges cause vibrations and large pressure forces, which in turn lead to broken preparation edges and rough surfaces. Instruments that are bent and/or do not work properly should be discarded immediately. Damaged, corroded or worn devices must not come into contact with intact instruments to avoid contact corrosion.

OPERATING INSTRUCTIONS

Follow the steps below for the drilling process of the bits.

- Fit the Drill to the Contra-Angle and configure the surgical motor with drilling speed, as indicated by the selected implant.
- According to planning, take the drill to the location of the surgical socket.
- Start the motor and drill with continuous insertion and removal movements, with abundant irrigation. This irrigation can be either manual or combined with engine irrigation.
- During drilling, pressure must not be excessive.
- Drilling must occur according to the length of the implant and the laser marking of the drill.
- The indication of drill sequence and rotation speed for each implant must be respected, contributing to the

success of the implant's osseointegration.

- Do not stop the motor rotation with the drill inside the surgical cavity, as this may make removal difficult or cause the drill to fracture. During piercing, use lip balm.

MANUAL CLEANING AND DISINFECTION SANITATION

This product must be properly sanitized after each use. Proceed as follows:

1. Dismantle the instruments (when applicable).
2. Immerse the instruments for at least 1 minute in the enzymatic detergent (CIDEZYME ®, 1.6% v/v) so that the instruments are sufficiently covered. Be careful that there is no contact between the instruments.
3. Carefully use a soft brush to help with the cleaning process. Shake instruments several times during cleaning.
4. Immerse the instruments for 15 minutes in the cleaning solution (CIDEZYME ®, 1.6% v/v) under ultrasonic treatment, so that the instruments are sufficiently covered. Be careful that there is no contact between the instruments.
5. Remove the instruments from the cleaning solution and wash them well for at least 3 times (for at least 1 minute) under running water.

DISINFECTION

1. Immerse the instruments (disassembled, if applicable) for 10 minutes in the disinfectant solution (CIDEX® OPA - OPA Solution - undiluted) so that the instruments are sufficiently covered.
2. Remove the instruments from the disinfectant solution and wash them according to the instructions below:

WASHING INSTRUCTIONS

1. After removing the instruments from the CIDEX® OPA Solution - OPA Solution, wash the medical device thoroughly by immersing it completely in a large volume of water. Use sterile water.
2. Keep the device fully immersed for at least 1 minute.
3. Remove the device and discard the wash water.
4. Always use new volumes of water for each wash.
5. Repeat the procedure 2 more times, for a total of 3 WASHES, with large volumes of clean water to remove residues of the CIDEX® OPA SOLUTION - OPA (Residues can cause serious side effects)
6. Inspect and package instruments immediately after removal.

AUTOMATIC CLEANING AND DISINFECTION

1. Use Neodisher® MediZym detergent.
2. Dismantle the instruments, if necessary;
3. Transfer the instruments to the Disinfectant Washer (be careful that the instruments do not come into contact);
4. Start the program;
5. Remove the instruments from the Disinfectant Washer after the end of the program;
6. Check and package instruments immediately after removal.

STERILIZATION

This product is reusable, supplied non-sterile and packaged individually.

This product must be properly sanitized and sterilized before use.

Sterilize it the day before or on the day of the procedure.

ATTENTION: These products must not be autoclaved in their original packaging.

For sterilization, only use the steam sterilization method according to the parameters below:

	Fractional Vacuum/Dynamic air Removal ¹	Gravitational ²
Sterilization Time	4 minutes	15 minutes
Sterilization Temperature³	132°C / 270°F	132°C / 270°F
Drying Time	At least 20 minutes ⁴	At least 20 minutes ⁴

1. At least three vacuum steps.
2. The less effective gravitational sterilization procedure should not be used if the fractional vacuum procedure is available.
3. Maximum sterilization temperature 134°C (273°F). The required effectiveness in drying time depends directly on the parameters under the user's responsibility (configuration and load density, sterilization conditions, and these must be determined by the user. However, the applied drying time must not be less than 20 minutes).

NOTES:

1. After sterilization, pack the instruments in a dry, dust-free environment.
2. The immediate/flash sterilization procedure should not be used.
3. Do not use thermal dry sterilization, radiation sterilization, formaldehyde and ethylene oxide sterilization, as well as plasma sterilization.

PRECAUTIONS

The Surgical Instrument Kit is supplied in non-sterile packaging. It is up to the team to sterilize the product before use, following classic autoclaving and biosafety protocols.

ADVERSE EFFECTS

No adverse effects are expected, as long as the product is used in accordance with the instructions for use.

ADDITIONAL INFORMATION FOR THE PROFESSIONAL

Instruct the patient about the need for professional medical monitoring after surgery and follow instructions on precautions, hygiene and medication prescription. These guidelines are the responsibility of the responsible professional.

USEFUL LIFE

This product is recommended for up to 20 uses, as long as the conditions of use recommended by DSP Biomedical are respected. Regardless of the number of times the instrument is used, the professional must always evaluate its conditions after each use.

STORAGE CONDITIONS

This product must be stored in its original packaging in a clean, dry place, with a maximum temperature of 45°C and protected from direct sunlight.

DISPOSAL OF MATERIAL

Any product and consumables used during surgery to install dental implants can compromise the health of those who handle them after use. Before disposing of them in the environment, it is recommended to observe current legislation and adhere to it.

GUARANTEE

DSP Biomedical guarantees the owner of this product a guarantee against any material or manufacturing defect. The presence of any defect must be communicated immediately to the manufacturer, respecting the legal deadline. The warranty for products manufactured by DSP Biomedical is linked to strict adherence to the information described in the instructions for use. Inappropriate use of the product, ignoring the instructions, exempts the manufacturer and/or distributor from any responsibility. Note: the warranty does not cover wear and tear of the product.





SUPPORT INFORMATION











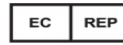

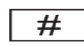







If there is a need for further information, or if the product presents any adverse effect, with potential risk to the patient, which generates or has the potential for injury or threat to public health, or any customer dissatisfaction, you must contact DSP through telephone 0800 600 88 66, or send an email to sac@dspbiomedical.com.br.

VALIDITY

8 years in sealed packaging unused

SÍMBOLOS

SIMBOLOGIA	DESCRIÇÃO	SIMBOLOGIA	DESCRIÇÃO
	Batch Number		See instructions for use
	Date manufacturing		Attention

	Manufacturing		Keep dry
	Non Sterile		Keep away from sunlight
	Product Code		Vality Date
	Temperature limit		Unique device identifier
	Umity limit		Do not use if the packaging is damaged and consult the instructions for use
	European Representative		Country of manufacturer
	Model Number		Medical Device
	Fragile, handle with care		Importer
	CE Mark		CE marking with number of Notified body; SIQ, number 1304
	Do not reuse		FDA-mandated prescription notification for the United States market

FABRICADO POR

DSP INDUSTRIAL LTDA
Rua Marechal Floriano Peixoto, 303 – Ouro Verde II
Campo Largo /PR – Brasil
CNPJ 03.960.018/0001-23
Telefone: +55 41 3291-2200
www.dspbiomedical.com
Responsável Técnico: CREA- PR 25412/D

REPRESENTANTE NA COMUNIDADE EUROPEIA

DSP BIOMEDICAL EUROPA UNIP LDA
Alameda dos Oceanos, 142 Lt. 4.24 0H
Parque das Nações – Lisboa - Portugal
1990-502
Telefone: (351) 962833592

Registro da ANVISA: 80116980004

*Federal (USA) law restricts this device to sale by or on the order of a licensed dentist or physician. Not all products are available in all countries. Please, contact the authorized distributor.